

Scheme 2026, Assessment Run 78

Aalborg, 2nd September 2026

Dear participant

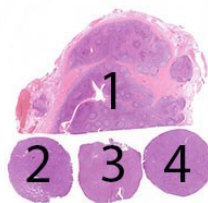
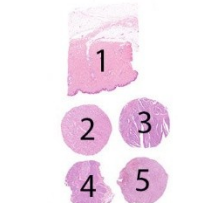
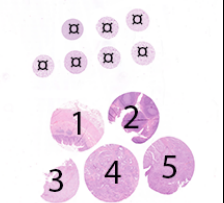
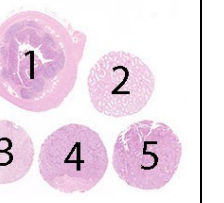
Pairs of unstained slides of NordiQC multi-tissue samples have been sent according to your submitted protocols for this run.

If you are in any doubt, please read the details on the website: www.nordiqc.org, or send us an e-mail: nordiqc@rn.dk. Also inform us if slides are missing.

Please note your participant number (NQC...) shown on the address label. This number **must** be indicated on all slides and used for any communication with NordiQC.

Contents of the NordiQC multi-tissue samples

Top row closest to the glass label. Please note the shape of the cores may vary from block to block.

				
CD23: 1: Tonsil 2: Mantle cell lymphoma 3-4: B-CLL	INSM1: 1: Appendix 2: Pancreas 3: Lung adenocarcinoma 4: Small cell lung carcinoma 5: Neuroendocrine tumour	MLA: 1: Skin 2: Kidney 3: Colon adenocarcinoma 4-5: Melanoma	P53: □: Cell lines* 1: Appendix 2: Tonsil 3: Low grade endometrial carcinoma 4-5: Endometrial serous carcinomas	SALL4: 1: Appendix 2: Testis 3: Renal clear cell carcinoma 4: Seminoma 5: Embryonal carcinoma

**Please note that the assessment slides for p53 contain additional cell lines beyond the NordiQC tissue cores. These are for research purposes only and are not part of the NordiQC assessment.*

Data generated from these additional cell lines may be used in future research publications.

For each marker stain the sections according to the protocols you have entered on our website.

If any changes are made in the protocols, update the used protocols as these are used for data analysis and individual feedback in case needed.

In-house controls shall, if possible, be placed on the NordiQC slide (see next page).

Mark all slides with your NordiQC participant number.

For each marker return only one stained slide. Keep the remaining NordiQC slides to serve as own control set. Do **not** return in-house controls on separate slides.

Deadline for slide return: **10th October.**

Slides return address:

NordiQC

Institute of Pathology

Aalborg University Hospital

Hospitalsbyen 22, Port 99.20 – 99.26

DK-9260 Gistrup

Denmark

Overall and individual results are available on the website on 17th December.

Thank you for your assistance.

Best regards
NordiQC

General information;

NordiQC tissue material and sections:

All included tissues in the NordiQC multi-tissue samples have been fixed in 10% neutral buffered formalin and embedded in paraffin.

The multi-tissue blocks are sectioned at 3-4 µm.

After sectioning, the slides are air dried at room temp. for 48 hours, after which they are stored at -80°C until distributed.

After receipt of slides, store them at -20°C until staining.

The slides are unbaked.

Before IHC staining in-house control material should be mounted on the NordiQC slide and dried at e.g. 37°C overnight or 55-60°C for one hour to secure section adhesion.

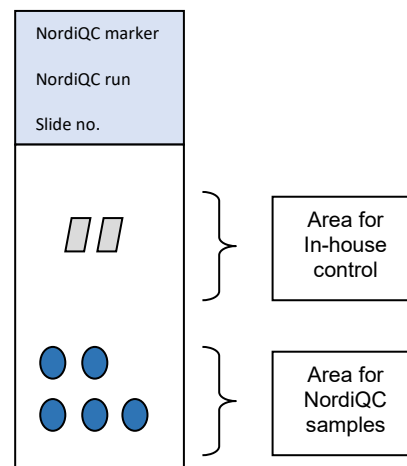


Fig. 1 Example of NordiQC slide

Slide return

Return slides to NordiQC after IHC staining has been completed and in due time for deadlines indicated in this accompany letter and at www.nordiqc.org/modules.php

The slides must be returned in the same or a similar box in which it was received.

You can see if NordiQC has received your slides. Please log in at our website and choose "Protocol submission". When the slides are received at NordiQC and registered, this will be marked with the date received in the column "Slide received by NordiQC".

Submitted sections

Submitted sections are stored by NordiQC for future documentation and will not be returned to the laboratories. However, laboratories can request their stained sections for review e.g. in case their own control set is lost.

In-house on-slide controls

In-house on-slide controls shall serve as participants internal IHC system control to verify successful test.

In-house on-slide controls will not be evaluated by NordiQC as expected patterns and content of tissue are not known and validated by NordiQC. The results in the controls will not affect the final individual assessment mark.

Results

The overall assessment results for all markers are published in reports available at www.nordiqc.org/epitope at the time indicated at <https://nordiqc.org/modules.php>.

The assessment reports will provide an overview of assessment criteria, protocols, reagents and instruments used by the participants and results obtained. The report will provide guidelines for recommendable control materials and descriptions / illustrations of the expected staining patterns in these.

Protocols giving optimal results will be published on the NordiQC website. Best practice protocols will typically and if possible be listed for the main IHC systems used by the participants.

In case of an insufficient result assessed as borderline or poor, the participant will receive individual suggestions for protocol improvement.