

## Scheme 2022, Assessment Run H22

Aalborg, 5<sup>th</sup> September 2022

Dear participant


Pairs of unstained slides of NordiQC multi-tissue samples have been sent according to your submitted protocols for this run.

If you are in any doubt, please read the details on the website: [www.nordiqc.org](http://www.nordiqc.org), or send us an e-mail: [nordiqc@rn.dk](mailto:nordiqc@rn.dk). Also inform us if slides are missing.

Please note your participant number (NQC...) shown on the address label. This number **must** be indicated on all slides and used for any communication with NordiQC.

### Contents of the NordiQC multi-tissue samples

Top row closest to the glass label.

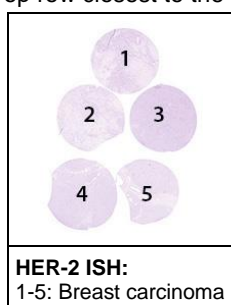

**IMPORTANT INFORMATION: Due to change of slide type, the writing on the slides can disappear during HIER. Use of a label can prevent this.**

Stain the section according to the protocol you have uploaded on our website. If any changes are made in the protocol, upload the used protocol.

In-house controls shall, if possible, be placed on the NordiQC slide (see next page).

Score the slide using the electronic scoring sheet on the web site (login and choose "scoring sheet" in the "options"-menu on the right).

HER-2 FISH: Do not return the FISH-stained slide to NordiQC.

HER-2 BRISH: Mark the slides with your NordiQC participant number and return one stained slide to NordiQC. Keep the remaining NordiQC slide to serve as a control set.

Do not return separate in-house controls.

Deadline for slide return: **10<sup>th</sup> October**

Slides return address:  
NordiQC  
Institute of Pathology  
Aalborg University Hospital  
Ladegaardsgade 3  
9000 Aalborg  
Denmark

Overall and individual results are available on the website on 17<sup>th</sup> December.

Thank you for your assistance.

Best regards  
NordiQC

General information;

### NordiQC tissue material and sections:

All included tissues in the NordiQC multi-tissue samples have been fixed in 10% neutral buffered formalin and embedded in paraffin.

The sections of the multi-tissue blocks are sectioned at 3-4  $\mu\text{m}$ .

After sectioning, the slides are air dried at room temp. for 48 hours, after which they are stored at  $-80^{\circ}\text{C}$  until distributed.

After receipt of slides, store them at  $-20^{\circ}\text{C}$  until staining.

The slides are unbaked.

Before IHC staining in-house control material should be mounted on the NordiQC slide and dried at e.g.  $37^{\circ}\text{C}$  overnight or  $55-60^{\circ}\text{C}$  for one hour to secure section adhesion.

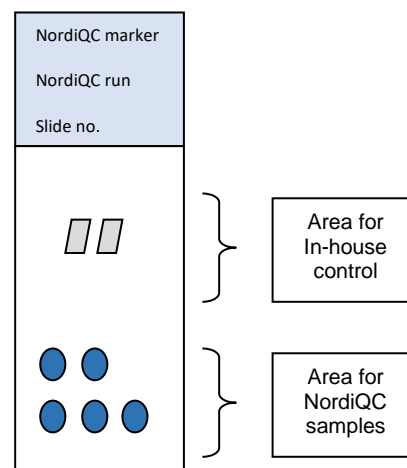


Fig. 1 Example of NordiQC slide

### Slide return

Return slides to NordiQC after BRISH staining has been completed and in due time for deadlines indicated in this accompany letter and at [www.nordiqc.org/modules](http://www.nordiqc.org/modules).

You can see if NordiQC has received your slides. Please log in at our website and choose "Protocol submission". When the slides are received at NordiQC and registered, this will be marked with the date received in the column "Slide received by NordiQC".

### Submitted sections

Submitted sections are stored by NordiQC for future documentation and will not be returned to the laboratories. However, laboratories can request their stained sections for review e.g. in case their own control set is lost.

### Results

The overall assessment results for all markers are published in reports available at [www.nordiqc.org/epitope](http://www.nordiqc.org/epitope) at the time indicated at <https://nordiqc.org/modules.php>.

The assessment reports will provide an overview of assessment criteria, protocols, reagents and instruments used by the participants and results obtained. The report will provide guidelines for recommendable control materials and descriptions / illustrations of the expected staining patterns in these.

In case of an insufficient result assessed as borderline or poor, the participant will receive individual suggestions for protocol improvement.